

AD-A251 078



USAARL Report No. 92-20

DTIC

ELECTF

JUN 4 1992

S C D

①



**Test and Evaluation Report
of the
Ohio Medical Transport Incubator Model Air-Vac**

By

James E. Bruckart (Project Officer)

Martin Quattlebaum (Project Officer)

Joseph R. Licina (Project Officer)

Biodynamics Research Division

Bill Olding (UES, Inc.)

92-14609



92 6 03 021

April 1992

Approved for public release; distribution unlimited.

**United States Army Aeromedical Research Laboratory
Fort Rucker, Alabama 36362-5292**

Notice

Qualified requesters

Qualified requesters may obtain copies from the Defense Technical Information Center (DTIC), Cameron Station, Alexandria, Virginia 22314. Orders will be expedited if placed through the librarian or other person designated to request documents from DTIC.

Change of address

Organizations receiving reports from the U.S. Army Aeromedical Research Laboratory on automatic mailing lists should confirm correct address when corresponding about laboratory reports.

Disposition

Destroy this document when it is no longer needed. Do not return it to the originator.

Disclaimer


The views, opinions, and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy, or decision, unless so designated by other official documentation. Citation of trade names in this report does not constitute an official Department of the Army endorsement or approval of the use of such commercial items.

Reviewed:



ROBERT W. WEIEN
MAJ, MC, SFS
Director, Biodynamics
Research Division

Released for publication:

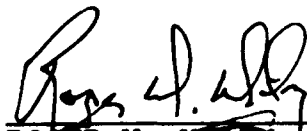

ROGER W. WILEY, O.D., Ph.D.
Chairman, Scientific
Review Committee



DAVID H. KARNEY
Colonel, MC, SRS
Commanding

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE

REPORT DOCUMENTATION PAGE				Form Approved OMB No. 0704-0188	
1a. REPORT SECURITY CLASSIFICATION UNCLASSIFIED			1b. RESTRICTIVE MARKINGS		
2a. SECURITY CLASSIFICATION AUTHORITY			3. DISTRIBUTION/AVAILABILITY OF REPORT Approved for public release, distribution unlimited		
2b. DECLASSIFICATION/DOWNGRADING SCHEDULE					
4. PERFORMING ORGANIZATION REPORT NUMBER(S) USAARL Report No. 92-20			5. MONITORING ORGANIZATION REPORT NUMBER(S)		
6a. NAME OF PERFORMING ORGANIZATION U.S. Army Aeromedical Research Laboratory		6b. OFFICE SYMBOL (If applicable) SGRD-UAD-IE	7a. NAME OF MONITORING ORGANIZATION U.S. Army Medical Research and Development Command		
6c. ADDRESS (City, State, and ZIP Code) P.O. Box 577 Fort Rucker, AL 36362-5292			7b. ADDRESS (City, State, and ZIP Code) Fort Detrick Frederick, MD 21702-5012		
8a. NAME OF FUNDING/SPONSORING ORGANIZATION		8b. OFFICE SYMBOL (If applicable)	9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER Partial effort under contract No. DAMD 17-86-C-6215		
8c. ADDRESS (City, State, and ZIP Code)			10. SOURCE OF FUNDING NUMBERS		
PROGRAM ELEMENT NO. 0603807A		PROJECT NO. 3M463807D8	TASK NO. 36 LC	WORK UNIT ACCESSION NO. 201	
11. TITLE (Include Security Classification) Test and Evaluation Report of the Ohio Medical Transport Incubator Model Air-Vac					
12. PERSONAL AUTHOR(S) James E. Bruckart, Martin Quattlebaum, Joseph R. Licina, and Bill Olding					
13a. TYPE OF REPORT Final		13b. TIME COVERED FROM _____ TO _____		14. DATE OF REPORT (Year, Month, Day) 1992 April	
15. PAGE COUNT 68					
16. SUPPLEMENTARY NOTATION					
17. COSATI CODES			18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number)		
FIELD	GROUP	SUB-GROUP	Electromagnetic compatibility, test and evaluation, aeromedical equipment		
01	03				
14	02				
19. ABSTRACT (Continue on reverse if necessary and identify by block number) The Ohio Medical Transport Incubator, Model Air-Vac, was tested for electromagnetic interference/compatibility in the UH-60A helicopter under the U.S. Army Program for Testing and Evaluation of Equipment for Aeromedical Operations. The tests were conducted using current military and industrial standards and procedures for electromagnetic interference/compatibility and human factors. The Ohio Medical Transport Incubator, Model Air-Vac, was found to be compatible with U.S. Army MEDEVAC UH-60A Black Hawk.					
20. DISTRIBUTION/AVAILABILITY OF ABSTRACT <input checked="" type="checkbox"/> UNCLASSIFIED/UNLIMITED <input type="checkbox"/> SAME AS RPT. <input type="checkbox"/> DTIC USERS			21. ABSTRACT SECURITY CLASSIFICATION UNCLASSIFIED		
22a. NAME OF RESPONSIBLE INDIVIDUAL Chief, Scientific Information Center			22b. TELEPHONE (Include Area Code) (205) 255-6907		22c. OFFICE SYMBOL SGRD-UAX-SI

Table of contents

SECTION	PAGE
1. <u>EXECUTIVE DIGEST</u>	
1.1 Test objectives	1-1
1.2 Testing authority	1-2
1.3 Scope	1-2
1.4 Material description	1-3
1.5 Summary	1-3
1.6 Conclusion	1-5
2. <u>SUBTESTS</u>	
2.1 Initial inspection	2-1
2.2 Battery life evaluation	2-1
2.3 Electrical safety evaluation	2-2
2.4 Human factors evaluation (laboratory)	2-3
2.5 Altitude (low pressure) test	2-3
2.6 Vibration test	2-4
2.7 High temperature test	2-6
2.8 Low temperature test	2-7
2.9 Humidity test	2-8
2.10 Electromagnetic characteristics test	2-9
2.11 In-flight human factors evaluation	2-11
2.12 In-flight EMI/EMC characteristics test	2-12
3. <u>SUPPORTING DOCUMENTATION</u>	
3.1 Detailed test information	3-1
3.2 Test data	3-3

3.3	Criteria, significant problems, and	3-36
	suggested improvements	
3.4	References	3-39
3.5	Abbreviations	3-40
3.6	List of manufacturers.	3-42
3.7	Distribution list	3-43

Section 1. Executive digest

The Army program for Test and Evaluation of Aeromedical Equipment uses existing military standards (MIL-STD) and collective professional expertise to test and evaluate selected medical equipment proposed for use aboard Army aircraft. Equipment meeting these standards ensures the safety of the crew, patients, and aircraft by eliminating risks due to: (1) Interference by the medical equipment with aircraft systems/subsystems operation, (2) the aircraft system's interference with the operation of the medical equipment, (3) the medical equipment's susceptibility to environmental exposure, or (4) physical and/or functional incompatibility while in use on board selected rotary-wing aircraft. This program tests both developmental and nondevelopmental (off the shelf) medical equipment destined for use aboard Army medical evacuation (MEDEVAC) aircraft.

1.1 TEST OBJECTIVES

1.1.1 To determine if the medical equipment is complete and operational per the manufacturer's operating instructions.

1.1.2 To ensure the electrical safety of the medical equipment.

1.1.3 To ensure the equipment will function as designed throughout the rated battery operation time.

1.1.4 To ensure the safety of the operator, the patient, and the aircrew.

1.1.5 To assess design considerations which could potentially contribute to an operator error.

1.1.6 To determine if the medical equipment can function as designed in a low pressure environment.

1.1.7 To determine the ability of the medical equipment to withstand the vibrational stresses expected in a rotary-wing flight environment without degradation or malfunction.

1.1.8 To determine the ability of the medical equipment to be stored and operated in a high temperature environment.

1.1.9 To determine the ability of the medical equipment to be stored and operated in a low temperature environment.

1.1.10 To determine the ability of the medical equipment to operate satisfactorily for short periods during exposure to high humidity conditions.

Accession For
NTIS GRA&I
DTIC TAB
Unannounced
Justification

Distribution/
Availability

Avail and
Dist Special



A-1

1.1.11 To assess the levels of electromagnetic emissions produced by the medical equipment within selected frequency ranges.

1.1.12 To assess the minimum electromagnetic susceptibility levels of the medical equipment within selected frequency ranges.

1.1.13 To assess the physical and/or functional compatibility of the medical equipment while in use on board the aircraft.

1.1.14 To assess the electromagnetic interference (EMI) and electromagnetic compatibility (EMC) characteristics of the medical equipment with the host aircraft and its installed systems.

1.2 TESTING AUTHORITY

Research and Technology Work Unit Summary, dated 5 October 1989. Project number 3M463807D836, titled, Army Program for Testing and Evaluation of Equipment for Aeromedical Operations.

1.3 SCOPE

1.3.1 This test was conducted at the United States Army Aeromedical Research Laboratory (USAARL), Cairns Army Airfield (CAAF), and designated test flight areas in and around Fort Rucker, Alabama.

1.3.2 The USAARL UH-60A aircraft, serial number 88-26069, with subsystems delineated in paragraph 3.2.2, was configured with the Ohio Medical Transport Incubator*, model Air-Vac and used as the test aircraft for the in-flight evaluation. The in-flight evaluation required 3.3 flight hours.

1.3.3 Laboratory testing was accomplished at USAARL using government furnished equipment (GFE) by Universal Energy Systems, Inc. (UES), under contract No. DAMD 17-86-C-6215.

1.3.4 Prior to flight testing, the following tests were accomplished: Acceptance inspection, equipment training, electromagnetic compatibility, human factors and safety, environmental compatibility, and in-flight compatibility.

1.3.5 An airworthiness release (AWR) dated 24 Feb 1992 was received from the U.S. Army Aviation Systems Command (AVSCOM) prior to the in-flight testing of the Ohio Transport Incubator.

* See list of manufacturers

1.4 MATERIAL DESCRIPTION

The Ohio Medical Air-Vac is an incubator system designed for intensive care isolation of infants during transportation on its mobile stand, in an ambulance, or in an aircraft. The incubator may be operated from 12 and 24 volts dc and 120 volts ac, 50 - 400 Hz power sources.

The incubator has a hinged plexiglass hood which covers the top and front of the infant compartment. A thermometer holder is located on the inside of the hood to measure internal temperature. Two 6-inch round access doors are located on the front of the hood to allow access to the infant from the side. A 6 x 9 inch access door is located at the end of the incubator to allow access to the infant's head and shoulders. The incubator is equipped with a support board and velcro straps to secure the infant and an adjustable IV stand stored on the back of the incubator.

A heating coil warms the air in the incubator. The temperature is selected by a temperature adjustment dial. Indicator lamps are provided to display POWER or HEATER operation. When the safety limit temperature is exceeded, a HIGH TEMP light is illuminated. Humidity may be added to the heated air via a humidity chamber on the side of the unit. Water is added to a sponge in a drawer and heated air flowing past the sponge provides 40-60% humidity for 3-4 hours. Provision is made for a "D" or "E" size oxygen cylinder in the body of the incubator to provide a portable source of oxygen. Oxygen concentrations of 40-90% may be obtained using an accessory regulator and fresh air intake vane on the incubator.

1.5 SUMMARY

1.5.1 Laboratory testing

1.5.1.1 Battery Life Evaluation: The Ohio Incubator used for testing was not factory new. The used battery would not hold a charge and was not tested. The manufacturer specifies a battery life up to 3 hours while maintaining 90°F in an ambient temperature of 70°F or 1 hour in ambient temperature of 40°F.

1.5.1.2 Electrical Safety Evaluation: All measurements were within acceptable limits. No unsafe qualities were found in the Ohio Medical Air-Vac. The limits for currents and resistances were in accordance with (IAW) the limits specified in TB-38-750-2, April 1987 and National Fire Prevention Association (NFPA) standards.

1.5.1.3 Human Factors Evaluation: The Ohio Medical Air-Vac was found to be satisfactory in all categories of the evaluation.

1.5.1.4 Environmental Tests: The Ohio Medical Air-Vac can be expected to perform in a variety of environmental conditions. Its performance was found to be satisfactory in all stages of the environmental testing except the high temperature and low temperature operation tests. The unit tested in the laboratory was in use for several years. When placed in a high temperature, the high temperature alarm did not activate. This indicated the high temperature alarm limit on this unit was not calibrated correctly at 100°F. The unit was not able to maintain 90°F inside the incubator when the ambient temperature was 32°F. The incubator temperature dropped to 73°F. The requirements for environmental tests are established in MIL-STD-810D, Methods 500.2 (altitude), 514.3 (vibration), 501.2 (high temperature), 502.2 (low temperature), and 507.2 (humidity).

1.5.1.5 Radiated Emissions Tests (RE02): The Ohio Medical Air-Vac may be unsatisfactory for use in certain EMI sensitive environments. Broadband (BB) radiated emissions were detected in the test frequency ranges. Some broadband emissions exceeded the test limits. Emission limits are set forth in MIL-STD-461A, Notice 4.

1.5.1.6 Radiated Susceptibility Test (RS03): The Ohio Medical Air-Vac was not found to be susceptible to radio frequency interference in the testing range and magnitude.

1.5.1.7 Conducted Emissions Test (CE01, CE02, and CE04): No signal failures were detected from the Ohio Transport Incubator during this test.

1.5.1.8 Conducted Susceptibility Test (CS02 and CS06): No susceptibility to the test power line spikes was noted in the Ohio Incubator.

1.5.2 In-flight testing

1.5.2.1 During the in-flight human factors evaluation, the Ohio Medical Air-Vac was found to be satisfactory in all categories of the evaluation criteria. The battery provided power for 45 minutes with the light on and the heater set for 90°F (58°F ambient temperature). If ac power is not available, the battery power may not be sufficient for a long transport in cold weather. The unit functioned properly in the aircraft with ac power.

1.5.2.2 The aircraft and its subsystems were not adversely affected by the operation of the Ohio Medical Air-Vac in any of the prescribed flight test modes.

1.5.2.3 The Ohio Medical Air-Vac was not affected by the aircraft and its subsystems during the in-flight testing.

1.6 CONCLUSIONS

Based on the results of laboratory and in-flight testing, the Ohio Medical Transport Incubator, Model Air-Vac was found to be compatible with U.S. Army MEDEVAC UH-60A Black Hawk with the subsystems listed in paragraph 3.2.2.

Section 2. Subtests

2.1 INITIAL INSPECTION

2.1.1 Objective

To determine if the Ohio Medical Air-Vac is complete and operational for testing per the manufacturer's operating instructions.

2.1.2 Criteria

2.1.2.1 The physical inventory is conducted solely for investigation and documentation.

2.1.2.2 The Ohio Medical Air-Vac will display consistent and accurate performance as an acceptable performance test.

2.1.3 Test procedure

2.1.3.1 A complete physical inventory of the Ohio Medical Air-Vac was completed per the manufacturer's equipment list.

2.1.3.2 An operational validation test of the Ohio Medical Air-Vac was conducted per the manufacturer's operating instructions by USAARL's medical maintenance personnel.

2.1.4 Test findings

2.1.4.1 The Ohio Medical Air-Vac was inventoried and found to be complete. The unit has been in service for several years prior to testing.

2.1.4.2 The Ohio Medical Air-Vac operated as prescribed in the manufacturer's operating manual. The used battery was not able to be recharged and a new replacement was not available. Criteria partially met.

2.2 BATTERY LIFE EVALUATION (Laboratory)

2.2.1 Objective

To ensure the equipment will function as designed throughout the rated battery operation time.

2.2.2 Criterion

Verify manufacturer's specified full power internal battery life expectancy of 3 hours operation while maintaining 90°F

internal temperature at 70°F ambient temperature or 1 hour operation in 40°F ambient temperature.

2.2.3 Test procedure

2.2.3.1 Charging and operation cycles were conducted in ambient room conditions.

2.2.4 Test findings

The battery indicated discharged after a full charging cycle so no test of battery life could be completed. The battery on a second unit was evaluated and provided 45 minutes service with light on and heater at 90°F (58°F ambient). The unit was not available for observation during 3 charge and discharge cycles. Criterion not evaluated.

2.3 ELECTRICAL SAFETY EVALUATION

2.3.1 Objective

To ensure the electrical safety, by evaluation of case-to-ground resistance and case-to-ground current leakage, of the Ohio Medical Air-Vac.

2.3.2 Criterion

The Ohio Medical Air-Vac shall meet the standards established in TB-38-750-2 and NFPA 99 for electrical safety of medical equipment.

2.3.3 Test procedure

Performance in the electrical safety evaluation were made, with a Neurodyne-Dempsey model 431F electrical safety analyzer*, IAW the procedures described in Technical Bulletin (TB) Number 38-750-2. Case-to-ground resistance and various case-to-ground leakage currents were measured. Leakage currents were measured using a 10 by 20 centimeter (cm) aluminum foil sheet taped flush to the equipment case. Checks were made for safety concerns such as case integrity, breaks in power cord insulation, and connectors.

2.3.4 Test findings

Grounding conductor resistance was 76.7 milliohms and maximum case leakage current was 10.7 microamperes. These measurements are below the limits specified in TB-38-750-2 and NFPA 99. Some cracked insulation on the power cords was noted and corrected. This deficiency would not be expected on a new unit. Criterion met.

2.4 HUMAN FACTORS EVALUATION (Laboratory)

2.4.1 Objectives

2.4.1.1 To assure the safety of the operator, the potential patient, and the aircrew.

2.4.1.2 To assess the design considerations which could potentially contribute to an operator error.

2.4.2 Criterion

The Ohio Medical Air-Vac must be rated satisfactory in all major categories of the evaluation. These include visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.

2.4.3 Test procedure

2.4.3.1 The evaluation was conducted in a laboratory under fluorescent lighting and ambient room conditions.

2.4.3.2 The Ohio Medical Air-Vac was operated according to prescribed instructions through its full range of functions.

2.4.4 Test finding

The first Ohio Medical Air-Vac was found to be satisfactory in all of the evaluation criteria except Conductors, Fasteners, and Safety. The power cord insulation was dry and cracked and the gas transport tubes are deteriorated. Some access door latches were broken and screws missing from the handle assemble. The operator's manual describes a battery pack that is different than the one supplied with this unit. These deficiencies would not be expected in a new unit. A second (newer) unit was evaluated and these deficiencies were not present. Criterion met.

2.5 ALTITUDE (LOW PRESSURE) TEST [IAW MIL-STD-810D, METHOD 500.2]

2.5.1 Objective

To determine if the Ohio Medical Air-Vac can function as designed in a low pressure environment.

2.5.2 Criterion

The Ohio Medical Air-Vac will perform as designed while exposed to an altitude equivalency of 15,000 feet above sea level.

2.5.3 Test procedure

2.5.3.1 A pretest performance check was conducted to ensure proper operation of the Ohio Medical Air-Vac.

2.5.3.2 The altitude test was performed in a Tenney Engineering model 64S altitude chamber*. This test is based on MIL-STD-810D, Method 500.2. The Ohio Medical Air-Vac was operated on the floor of the chamber. Chamber pressure was decreased to 420 mmHg (15,000 ft equivalent altitude) over a 15-minute period, held constant for 60 minutes, then raised, at 1500 fpm, to ambient conditions (760 mmHg) over a 10-minute period. There are no provisions for the control of temperature or humidity inside this chamber.

2.5.3.3 A posttest performance check was conducted to ensure proper operation of the Ohio Medical Air-Vac after the exposure to low pressure.

2.5.4 Test findings

2.5.4.1 The pretest performance check met criterion 2.1.2.2.

2.5.4.2 No failures in the performance of the Ohio Medical Air-Vac were noted before, during, or after the altitude test. Criterion met.

2.5.4.3 The posttest performance check met criterion 2.1.2.2.

2.6 VIBRATION TEST [IAW MIL-STD-810D, METHOD 514.3]

2.6.1 Objective

To determine the ability of the Ohio Medical Air-Vac to withstand the vibrational stresses expected in a rotary-wing environment without degradation or malfunction.

2.6.2 Criterion

The Ohio Medical Air-Vac will remain operational and be able to display consistent and accurate performance while exposed to vibrational stresses.

2.6.3 Test procedure

2.6.3.1 A pretest performance check was conducted to ensure proper operation of the Ohio Medical Air-Vac.

2.6.3.2 The vibration test was performed using an Unholtz-Dickey model TA115-40/CSTA vibration test system*. It is a single-axis system with an electromagnetic driver unit. The test consisted

of sinusoidal vibrations superimposed on random vibrations over a frequency range of 500 Hz, as shown below. These vibrations are derived from performance taken on the floor under the copilot's seat in a UH-1 helicopter traveling at 120 knots. The reference spectrum breakpoints are from MIL-STD-810D, Method 514.3; reference spectrum levels are based on field performance with a conservatism factor of 1.5. Independent tests were conducted in the X, Y, and Z axes.

Z-axis

duration: 60 minutes
broadband intensity: $0.4506 G_{rms}$
random vibration: initial slope : 99.00 dB/oct
5 Hz level: $0.00006210 G_{sqr/Hz}$
100 Hz level: $0.0006210 G_{sqr/Hz}$
300 Hz level: $0.0006210 G_{sqr/Hz}$
500 Hz level: $0.00006210 G_{sqr/Hz}$
final slope: -99.00 dB/oct
sinusoidal vibration: $.5450 G_{pk}$ at 11.25 Hz
 $.1690 G_{pk}$ at 22.50 Hz
 $.1200 G_{pk}$ at 33.75 Hz
 $.0310 G_{pk}$ at 45.00 Hz
 $.0530 G_{pk}$ at 56.25 Hz

X and Y axes

duration: 60 minutes each
broadband intensity: $0.3099 G_{rms}$
random vibration: initial slope: 99.00 dB/oct
5 Hz level: $0.00002920 G_{sqr/Hz}$
100 Hz level: $0.0002920 G_{sqr/Hz}$
300 Hz level: $0.0002920 G_{sqr/Hz}$
500 Hz level: $0.00002920 G_{sqr/Hz}$
final slope: -99.00 dB/oct
sinusoidal vibration: $.3200 G_{pk}$ at 11.25 Hz
 $.0670 G_{pk}$ at 22.50 Hz
 $.0950 G_{pk}$ at 33.75 Hz
 $.0350 G_{pk}$ at 45.00 Hz
 $.0770 G_{pk}$ at 56.25 Hz

The Ohio Medical Air-Vac was strapped to the vibration table fixture, and its performance was evaluated before, during, and after exposure to vibration.

2.6.3.3 A posttest performance check was conducted to ensure proper operation of the Ohio Medical Air-Vac.

2.6.4 Test findings

2.6.4.1 The pretest performance check met criterion 2.1.2.2.

2.6.4.2 No failures in the performance of the Ohio Medical Air-Vac occurred before, during, or after exposure to vibration. Criterion met.

2.6.4.3 The posttest performance check met criterion 2.1.2.2.

2.7 HIGH TEMPERATURE TEST [IAW MIL-STD-810D, METHOD 501.2]

2.7.1 Objective

To determine the ability of the Ohio Medical Air-Vac to be stored and operated in a high temperature environment.

2.7.2 Criteria

2.7.2.1 The Ohio Medical Air-Vac will demonstrate consistent and accurate operation during the high temperature operation check.

2.7.2.2 The Ohio Medical Air-Vac will demonstrate consistent and accurate operation after the high temperature storage cycle.

2.7.3 Test procedure

2.7.3.1 A pretest performance check was conducted to ensure proper operation of the Ohio Medical Air-Vac.

2.7.3.2 The high temperature test was conducted in a Tenney Engineering model ZWUL-10107D walk-in controlled environment chamber*. This test is based on MIL-STD-810D, Method 501.2. For the high temperature operation test, the Ohio Medical Air-Vac was turned on and placed on the floor of the environmental chamber. The chamber temperature was raised to 49°C and the humidity was stabilized at a maximum of 20 percent relative humidity (RH) within 15 minutes. The environmental control system is capable of regulating temperature within $\pm 2^\circ\text{C}$ and humidity within ± 5 percent RH. Temperature and humidity were held constant for 2 hours. At 30-minute intervals, the chamber door was opened briefly to minimize the change in chamber conditions during performance checks. After the operational test, the Ohio Medical Air-Vac was allowed to return to ambient conditions over a 30-minute period.

2.7.3.3 A posttest performance check was conducted to ensure proper operation of the Ohio Medical Air-Vac.

2.7.3.4 The Ohio Medical Air-Vac was stored (not operated) at temperatures of 63°C for 1 hour, 71°C for 4 hours, then again at 63°C for 1 hour. The chamber and Ohio Medical Air-Vac then were returned to ambient conditions over a 30-minute period.

2.7.3.5 A poststorage performance check was conducted to ensure proper performance of the Ohio Medical Air-Vac.

2.7.4 Test findings

2.7.4.1 The pretest performance check met criterion 2.1.2.2.

2.7.4.2 No operational failures occurred during the high temperature test. The high temperature alarm failed to activate when the internal temperature of the incubator exceeded 100°F. This indicates that the high temperature alarm limit was not properly calibrated. Criterion partially met.

2.7.4.3 The posttest performance check met criterion 2.1.2.2.

2.7.4.4 The Ohio Medical Air-Vac functioned properly after the high temperature storage test. Criterion met.

2.8 LOW TEMPERATURE TEST [IAW MIL-STD-810D, METHOD 502.2]

2.8.1 Objective

To determine the ability of the Ohio Medical Air-Vac to be stored and operated in a low temperature environment.

2.8.2 Criteria

2.8.2.1 The Ohio Medical Air-Vac will demonstrate consistent and accurate operation during the low temperature operation check.

2.8.2.2 The Ohio Medical Air-Vac will demonstrate consistent and accurate operation after the low temperature storage cycle.

2.8.3 Test procedure

2.8.3.1 A pretest performance check was conducted to ensure proper operation of the Ohio Medical Air-Vac.

2.8.3.2 The Ohio Medical Air-Vac was placed on the floor of the environmental chamber and the temperature was lowered to 0°C within 25 minutes. The environmental control system is capable of regulating temperature within 2°C. Humidity cannot be controlled in the chamber at freezing temperatures. The temperature was held constant for 2 hours. The chamber door was opened briefly every 30 minutes to minimize the change in chamber conditions, and a performance check was conducted. The chamber temperature then was raised to ambient temperature within a 30-minute period.

2.8.3.3 A posttest performance check was conducted to ensure proper operation of the Ohio Medical Air-Vac.

2.8.3.4 The Ohio Medical Air-Vac was "stored" in a nonoperational mode. The Ohio Medical Air-Vac was placed on the floor of the environmental test chamber and the temperature was lowered to -46°C for 6 hours. The chamber then was raised to ambient temperature over a 30-minute period.

2.8.3.5 A poststorage performance check was conducted to ensure proper operation of the Ohio Medical Air-Vac.

2.8.4 Test findings

2.8.4.1 The pretest performance check met criterion 2.1.2.2.

2.8.4.2 No operational failures occurred during the low temperature test. The incubator could not maintain the set temperature (90°F). With an ambient temperature of 32°F, the unit maintained an internal temperature of 73°F. Criterion partially met.

2.8.4.3 The posttest performance check met criterion 2.1.2.2.

2.8.4.4 The Ohio Medical Air-Vac functioned properly after the low temperature storage test. Criterion met.

2.9 HUMIDITY TEST [IAW MIL-STD-810D, METHOD 507.2]

2.9.1 Objective

To determine the ability of the Ohio Medical Air-Vac to operate satisfactorily for short periods of time during exposure to highly humid conditions.

2.9.2 Criterion

The Ohio Medical Air-Vac will demonstrate consistent and accurate operation while exposed to a high humidity environment.

2.9.3 Test procedure

2.9.3.1 A pretest performance check was conducted to ensure the proper operation of the Ohio Medical Air-Vac.

2.9.3.2 The humidity test was conducted in a Tenney Engineering model ZWUL-10107D walk-in controlled environment chamber*. This test is based on MIL-STD-810D, Method 507.2. For the humidity test, the Ohio Medical Air-Vac was placed in operation on the floor of the environmental chamber. The chamber temperature was raised to a temperature of 30°C and a relative humidity of 95 percent within 25 minutes. Temperature and relative humidity were maintained for 4 hours. The environmental control system is capable of regulating temperature within $\pm 2^\circ\text{C}$ and humidity within ± 5 percent RH. At 45-minute intervals the performance of

the Ohio Medical Air-Vac was checked. The chamber door was opened briefly to minimize the change in chamber conditions. The chamber and the Ohio Medical Air-Vac were returned to ambient conditions before the posttest performance validation check was conducted.

2.9.3.3 A posttest performance check was conducted to ensure the proper operation of the Ohio Medical Air-Vac.

2.9.4 Test findings

2.9.4.1 The pretest performance check met criterion 2.1.2.2.

2.9.4.2 No failures were noted in the Ohio Medical Air-Vac performance checks conducted during the exposure to the high humidity environment. Criterion met.

2.9.4.3 The posttest performance check met criterion 2.1.2.2.

2.10 ELECTROMAGNETIC CHARACTERISTICS TEST [IAW MIL-STD-461A, Notice 4, AND MIL-STD-462, Notice 3]

2.10.1 Objectives

2.10.1.1 To assess the maximum levels of radiated electromagnetic emissions produced by the Ohio Medical Air-Vac in the 14 kHz to 12.4 GHz frequency range.

2.10.1.2 To assess the tolerances of radiated electromagnetic susceptibility of the Ohio Medical Air-Vac within the 10 kHz to 10 GHz electric field.

2.10.1.3 To assess the maximum levels of conducted electromagnetic emissions produced by the Ohio Medical Air-Vac in the 10 kHz to 50 MHz frequency ranges.

2.10.1.4 To assess the tolerances of conducted electromagnetic susceptibility of the Ohio Medical Air-Vac within the range of 50 kHz to 400 MHz and power spikes.

2.10.2 Criteria

2.10.2.1 The Ohio Medical Air-Vac will not produce emissions in excess of the limits set forth in MIL-STD-461A, Notice 4, paragraph 6.13.

2.10.2.2 The Ohio Medical Air-Vac will not malfunction when it is subjected to radiated emissions as specified in MIL-STD-461A, Notice 4, paragraph 6.20.

2.10.2.3 The Ohio Medical Air-Vac will not conduct emissions in excess of the limits set forth in MIL-STD-461A, Notice 4, paragraphs 6.1 and 6.2.

2.10.2.4 The Ohio Medical Air-Vac will not malfunction when it is subjected to conducted emissions as specified in MIL-STD-461A, Notice 4, paragraphs 6.7 and 6.10.

2.10.3 Test procedure

2.10.3.1 The radiated emissions test was performed according to MIL-STD-462, Notice 3, Method RE02. The Ohio Medical Air-Vac was positioned on a wooden test stand inside the EMI chamber, 1 meter away from the receiving antennas. The antennas were mounted for both vertical and horizontal polarities and connected to EMI receivers. While the Ohio Medical Air-Vac was operating, the frequency spectrum (14 kHz to 12.4 GHz) was scanned for emissions. The Ohio Medical Air-Vac was operated with ac power.

2.10.3.2 The radiated susceptibility test was performed according to MIL-STD-462, Notice 3, Method RS03. The Ohio Medical Air-Vac was positioned on a wooden test stand inside the EMI chamber 1 meter away from the transmitting antennas. The antennas were mounted for both vertical and horizontal polarities and connected to radio frequency (RF) transmitters. While the Ohio Medical Air-Vac was operating, it was monitored for faulty operation during exposures to fields of 1 V/m from 10 kHz to 2 MHz, and 5 V/m from 2 to 30 MHz, 10 V/m from 30 MHz to 2 GHz, and 5 V/m from 2 to 10 GHz. The Ohio Medical Air-Vac was operated with ac power.

2.10.3.3 The conducted emissions tests were performed according to MIL-STD-462, Notice 3, Methods CE02 and CE04. The IMED Air-Vac was placed on a grounded, copper-covered workbench. The top of the workbench was 1 meter from floor level, 1.37 meters long and 0.81 meters wide. Power was supplied via a pair of line impedance stabilization networks (LISN) and a test jig. The test jig is a wooden tray with two power receptacles and two slots to hold current probes in place around power supply conductors. While the Ohio Medical Air-Vac was operating, the frequency range (10 kHz to 50 MHz) was scanned for emissions conducted in the power cable from the Ohio Medical Air-Vac.

2.10.3.4 The conducted susceptibility spike test was performed on a chemical resistant counter top according to MIL-STD-462, Notice 3, Method CS06. Power was supplied via a customized metal connection box. The connection box has two power receptacles and four banana jacks on its front panel. Connections to the individual power lines were made in series through the banana jacks. Transient spikes of 100 volts, 10 microseconds were generated with a Solar Electronics model 8282-1 transient pulse generator*

and induced onto the power leads at the connection box banana jacks. The spikes were monitored with a Tektronix 2235 oscilloscope* connected to a power receptacle on the connection box. The Ohio Medical Air-Vac was plugged into the other receptacle on the connection box and placed in operation. It was observed visually for correct operation while it was subjected to the power line spikes.

2.10.3.5 The conducted susceptibility test was performed according to MIL-STD-462, Notice 3, Method CS02. The Ohio Air-Vac was placed on a grounded, copper-covered workbench. Radio frequency interference was induced on the power leads and measured at the Ohio Air-Vac power cable. The frequency of the interference was incremented over the 50 kHz to 400 MHz range while the Ohio Air-Vac was operated. It was observed visually for proper operation while it was subjected to the radio interference on the power leads. Each frequency was held for 15 seconds.

2.10.4 Test findings

2.10.4.1 During the radiated emissions test, emissions which exceeded specification limits of MIL-STD-461A, Notice 4, were detected. These included:

<u>Frequency range</u>	<u>Emission exceeding standard</u>
150 - 245 kHz	0.4 - 1.7 dB (BB)
4.7 - 17.97 MHz	0.8 - 9.5 dB (BB)
23.44 - 27.81 MHz	1.7 - 13.4 dB (BB)

Criterion partially met.

2.10.4.2 The Ohio Medical Air-Vac was not found to be susceptible to radio frequency interference in the testing range and magnitude. Criterion met.

2.10.4.3 No signal failures were detected from the Ohio Air-Vac during the conducted emissions test. Criterion met.

2.10.4.4 The Ohio Air-Vac was not susceptible to radio frequency interference (RFI) or test spikes during the conducted susceptibility tests. Criterion met.

2.11 IN-FLIGHT HUMAN FACTORS EVALUATION

2.11.1 Objective

To assess the physical and/or functional compatibility of the Ohio Medical Air-Vac while in use onboard the aircraft.

2.11.2 Criterion

The flight surgeon will be able to operate the Ohio Medical Air-Vac without physical or functional restrictions aboard the aircraft. Major areas of concern include: Proper operation, visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.

2.11.3 Test procedure

2.11.3.1 A human factors evaluation was performed IAW MIL-STD-1472D, AAMI Human factors engineering guidelines, and UL-544 to ensure the compatibility of the Ohio Medical Air-Vac and the in-flight environment. The flight surgeon conducted the test wearing a flight suit, flight gloves, and an SPH-4B flight helmet. An evaluation of the compatibility with the nuclear, biological, and chemical (NBC) protective equipment was not conducted. Due to restrictions of the AWR, testing was conducted during daylight hours only.

2.11.3.2 The Ohio Medical Air-Vac was placed on the floor of the aircraft and secured with cargo straps. The Ohio Medical Air-Vac was tested using ac and battery power in all flight scenarios required by the In-flight Test Operations Procedures (ITOP) (refer to section 3.2).

2.11.4 Test findings

During the in-flight human factors evaluation, the Ohio Medical Air-Vac was found to be satisfactory in all categories of the evaluation criteria. The battery provided 45 minutes operation with the light on and the internal temperature set to 90°F (58°F ambient). Battery power alone may not be sufficient for a long transport or cold climate. Criterion met.

2.12 IN-FLIGHT EMI/EMC CHARACTERISTICS

2.12.1 Objective

To assess the EMI/EMC characteristics of the Ohio Medical Air-Vac with the host aircraft and its installed systems.

2.12.2 Criteria

2.12.2.1 The Ohio Medical Air-Vac will not radiate EMI to disrupt or interfere with other equipment or systems aboard the aircraft.

2.12.2.2 The aircraft will not radiate EMI to disrupt or interfere with the Ohio Medical Air-Vac's operation.

2.12.3 Test procedure

A qualitative EMI/EMC assessment was performed with both the Ohio Medical Air-Vac and the aircraft operating as source and victim. The Ohio Medical Air-Vac and applicable aircraft instruments and systems were monitored for unusual operation, readings, surges, or power anomalies for each checklist item (see pages 3-4 through 3-7).

2.12.4 Test findings

2.12.4.1 There were no adverse instances of EMI/EMC noted with the Ohio Medical Air-Vac acting as either the source or victim. Criterion met.

2.12.4.2 There were no adverse instances of EMI/EMC noted with the aircraft acting as either the source or victim. Criterion met.

Section 3. Supporting documentation

3.1 DETAILED TEST INFORMATION

3.1.1 General information

3.1.1.1 Ohio Medical Air-Vac testing is not considered a major action significantly affecting the quality of the human environment and, therefore, qualifies for categorical exclusion A-28, appendix A, AR 200-1.

3.1.1.2 A safety pilot will be designated for each flight. Flight operations will be conducted IAW the aircraft operator's manual, appropriate aircrew training manuals, and test item technical data.

3.1.2 Material description

3.1.2.1 The Ohio Medical Air-Vac is an incubator system designed for intensive care isolation of infants during transportation on its mobile stand, in an ambulance, or in an aircraft. The incubator may be operated from 12 and 24 volts dc and 120 volts ac, 50 - 400 Hz power sources.

The incubator has a hinged plexiglass hood which covers the top and front of the infant compartment. A thermometer holder is located on the inside of the hood to measure internal temperature. Two 6-inch round access doors are located on the front of the hood to allow access to the infant from the side. A 6 x 9 inch access door is located at the end of the incubator to allow access to the infant's head and shoulders. The incubator is equipped with a support board and velcro straps to secure the infant and an adjustable IV stand stored on the back of the incubator.

A heating coil warms the air in the incubator. The temperature is selected by a temperature adjustment dial. Indicator lamps are provided to display POWER or HEATER operation. When the safety limit temperature is exceeded, a HIGH TEMP light is illuminated. Humidity may be added to the heated air via a humidity chamber on the side of the unit. Water is added to a sponge in a drawer and heated air flowing past the sponge provides 40-60% humidity for 3-4 hours. Provision is made for a "D" or "E" size oxygen cylinder in the body of the incubator to provide a portable source of oxygen. Oxygen concentrations of 40-90% may be obtained using an accessory regulator and fresh air intake vane on the incubator.

3.1.2.2 Dimensions: 96.5 x 49.5 x 38.1 cm (38 x 19.5 x 15 in).


3.1.2.3 Weight: Incubator - 39.9 kg (88 lbs)
Transport cart - 16.8 kg (37 lbs)

3.1.2.4 Heating capacity: 12 or 24 Vdc - capable of maintaining 90°F in 35°F ambient; 120 Vac - capable of raising temperature 20°F above a 70°F ambient temperature in less than 20 minutes.

3.1.2.5 Power requirements: 120 Vac, 60 Hz, 2 amps; 240 Vac, 50 Hz, 1 amp; 12 Vdc, 14 amps; or 24 Vdc, 7 amps.

3.2 TEST DATA

3.2.1 Photographic description

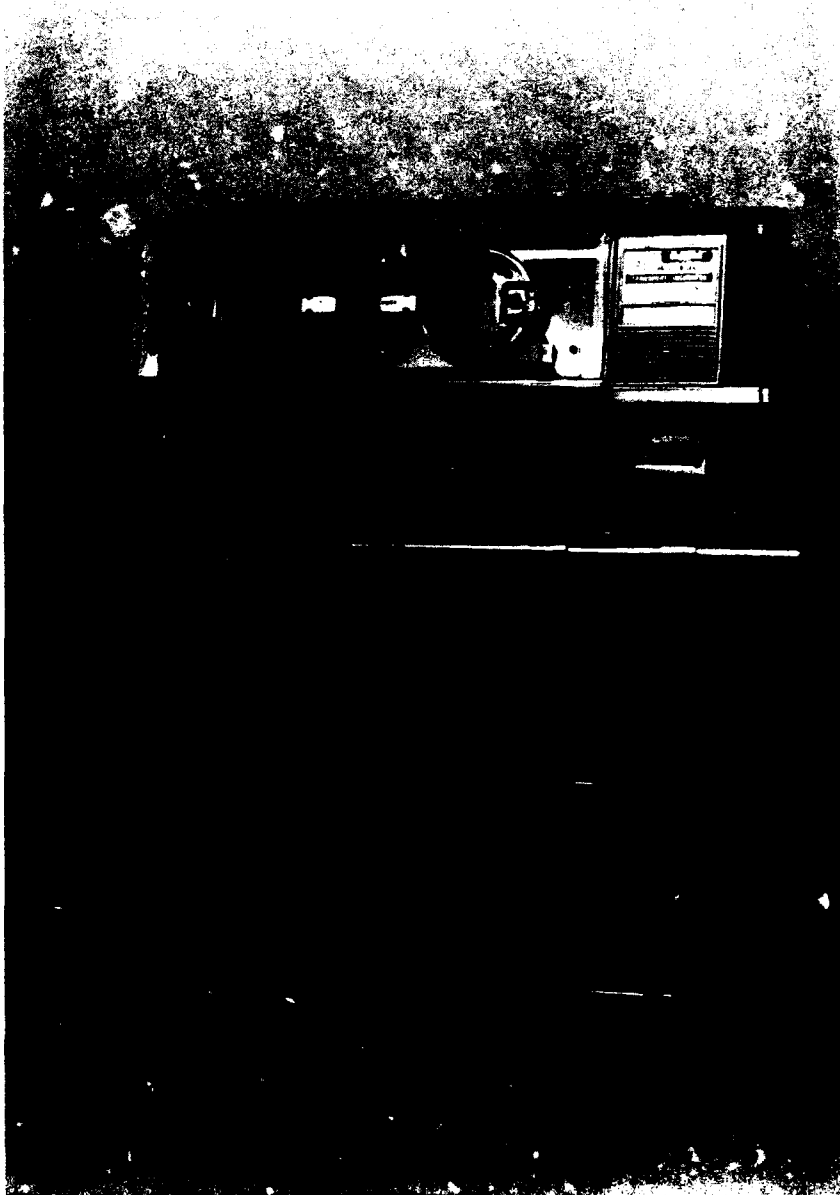

3.2.2 Aircraft equipment list

Item No.	Nomenclature
1	Receiver radio -- R-1496A/ARN-89 (automatic direction finder)
2	Displacement gyro -- CN-1314/A
3	Gyro directional -- CN-998/ASN-43
4	Signal data converter -- CV-3338/ASN-128
5	Receiver -- R-2139/ARN-123 (VOR/LOC/MB/GS)
6	Command instrument system processor -- 70600- 01038-101
7	SAS amplifier -- 70901-02908-104 (flight control stability augmentation system)
8	Rate gyro -- TRU-2A/A
9	Amplifier, impedance -- AM-4859A/ARN-89
10	Cargo hook -- FE-7590-145
11	Receiver, radar -- RT-1193/ASN-128 (doppler navigation receiver)
12	Barometric altimeter -- AAU-31/A-1
13	Barometric altimeter -- AAU-32A
14	Receiver/transmitter -- RT-1300/ARC-186 (VHF-AM and/or FM radio)
15	UHF-AM radio set -- RT-1518/ARC-164
16	Interphone control -- C6533/ARC (aircraft intercom control)
17	Receiver/transmitter -- RT-1115D/APN-209 (radar altimeter)
18	Indicator altimeter -- ID-1917C/APN-209 (radar altimeter)
19	Control radio set -- C-7392A/ARN-89 (automatic direction finder)
20	Comparator signal data -- CM-482/ARC-186 (comparator for ARC-186)
21	Receiver/transmitter -- RT-1296A/APX-100 (transponder with IFF)
22	Computer display unit -- CP-1252/ASN-128 (doppler navigation system)
23	Compass set controller -- C-8021E/ASN75
24	Magnetic compass - standby -- MS-17983-4

3.2.3 In-flight test data card

DATA CARD FORMAT

GUIDELINE FOR DATA COLLECTION

IN-FLIGHT SUITABILITY TEST OF MEDICAL ITEMS

1. Installation/removal.	Suitable		Comments
	Yes	No	
a. Weight and balance (DD Form 365-4, Clearance Form F).	X		
b. Space/area allocation.			
(1) Operational requirements.	X		
(2) Storage requirements.	X		
c. Interface connections (safe, positive, secure).	X		
d. Installation/removal (expedient/easily achieved).	X		
e. Mounting/final config- uration (functional/stable).	X		
2. Operations and performance.	Suitable		Comments
	Yes	No	
a. Manufacturer's operating instruction.	X		
b. Medical item operation before aircraft run-up.	X		
c. System interface during aircraft engine run-up and medical item operation (EMI switchology checklist).	X		
(1) Aircraft voltage output.	X		

	Suitable		Comments
	Yes	No	
(2) Flight control function (UH-60).	X		
(3) Stabilator function (UH-60).	X		
(4) Radio communication vs. medical item operation.			
(a) FM	X		
(b) UHF	X		
(c) VHF	X		
(5) Navigation equipment vs. medical item operation.			
(a) Transponder	X		
(b) ADF	X		
(c) VOR	X		
(d) Doppler	X		
(6) Radar altimeter operation vs. medical item operation.	X		
d. System interface during aircraft hover and medical item operation (EMI switchology checklist).			
(1) Voltage output.	NA		
(2) Radio communication vs. medical item operation.			
(a) FM	X		
(b) UHF	X		
(c) VHF	X		

(3) Navigation equipment operation vs. medical item operation.	Suitable		Comments
	Yes	No	
(a) Transponder	X		
(b) ADF	X		
(c) VOR	X		
(d) Doppler	X		
e. Flight mission profile vs. medical item operation (EMI switchology checklist).			
(1) Straight and level (1000 ft MSL for 20 minutes).			
(a) Compatibility of flight mode and medical item operation.	X		
(b) Radio communication vs. medical item opera- tion.			
<u>a.</u> FM	X		
<u>b.</u> UHF	X		
<u>c.</u> VHF	X		
(2) NOE (20 minutes). compatibility of flight mode and medical item operation.	X		
(3) FM homing (10 minutes).	X		
(4) Doppler navigation vs. medical item operation.			
(a) Initialize function.	X		
(b) Fix function.	X		
(c) Update function.	X		

	Suitable Yes No	Comments
(5) VOR navigation 7000 ft MSL for 20 minutes) vs. medical item operation.	X	
(6) ILS approach vs. medical item operation.	X	
f. Medical item operation after engine shutdown (external power source).	X	
g. Restrictions to the medical item's use (i.e., electrical connectors).	X	
h. Deviations from the labor- atory test results.		
(1) Electrical/ electronic.	None	
(2) Mechanical environment.	None	
(3) Human factors (user interface, controls, markings, lighting, egress).	None	
(4) Safety.	None	
3. Deviations from the in-flight test protocol.		
a. The VOR navigation portion of the in-flight test con- ducted at 2000 feet MSL due to air traffic control clearance.		

3.2.4 EMI switchology checklist

EMI SWITCHOLOGY CHECKLIST UH-60 AIRCRAFT

IN-FLIGHT SUITABILITY OF MEDICAL ITEMS

ENG INSTRUMENTS/CDU	No EMI Affect	EMI Affected Gnd Flt	Explanation
Fuel quantity	X		
Fuel indicator test	X		
XMSN oil temperature	X		
XMSN oil pressure	X		
#1 engine oil temperature	X		
#2 engine oil temperature	X		
#1 engine oil pressure	X		
#2 engine oil pressure	X		
#1 TGT	X		
#2 TGT	X		
#1 Ng speed	X		
#2 Ng speed	X		
CDU digits on/off	X		
CDU instruments dim	X		
ENG INSTRUMENTS/PLT PDU	No EMI Affect	EMI Affected Gnd Flt	Explanation
#1 engine RPM	X		
#2 engine RPM	X		
Rotor RPM	X		
#1 torque	X		
#2 torque	X		
ENG INSTRUMENTS/COPLT PDU	No EMI Affect	EMI Affected Gnd Flt	Explanation
#1 engine RPM	X		
#2 engine RPM	X		
Rotor RPM	X		
#1 torque	X		
#2 torque	X		

ENG CONTROLS	No EMI Affect	EMI Affected Gnd Flt	Explanation
#1 overspeed	X		
#2 overspeed	X		
RPM switch	X		
#1 engine anti-ice	X		
#2 engine anti-ice	X		
#1 inlet anti-ice	X		
#2 inlet anti-ice	X		

RADIO EQUIPMENT	No EMI Affect	EMI Affected Gnd Flt	Explanation
ICS, C-6533 ARC	X		
VHF-FM, ARC-186/115	X		
VHF-AM, ARC-186/115	X		
UHF-AM, ARC-164(V)	X		
Crypto, KY-28			Not installed
Radio retransmissions PLN			Not installed
Transponder, APX-100(V)	X		
KIT-1A/TSEC IFF computer			Not keyed with code

MISSION EQUIPMENT	No EMI Affect	EMI Affected Gnd Flt	Explanation
RWR, APR-39(V)			Not installed
IR CM, ALQ-144			Not installed
Chaff dispenser, M-130			Not installed
Cargo hook system	X		

HYDRAULIC CONTROL SYSTEM	No EMI Affect	EMI Affected Gnd Flt	Explanation
Backup hydraulic pump	X		
Servo off 1st stage/PLT	X		
Servo off 2nd stage/PLT	X		
Servo off 1st stage/COPLT	X		
Servo off 2nd stage/COPLT	X		
Hydraulic leak test	X		
Tail servo	X		
Boost servos	X		

FUEL SYSTEM	No EMI Affect	EMI Affected Gnd Flt	Explanation
Fuel pump switch	X		
Fuel boost pump #1	X		
Fuel boost pump #2	X		
Fuel cont panel ESSS	X		
WARNING SYSTEM	No EMI Affect	EMI Affected Gnd Flt	Explanation
Low rotor RPM	X		
Master caution	X		
Caution advisory	X		
Fire warning	X		
AFCS	X		
Stabilator	X		
#1 engine out	X		
#2 engine out	X		
NAVIGATION INSTRUMENTS	No EMI Affect	EMI Affected Gnd Flt	Explanation
ADF	X		
Magnetic compass	X		
CONUS NAV, ARN-123	X		
Doppler, ASN-128	X		
Gyro mag compass (PLT)	X		
Gyro mag compass (COPLT)	X		
Compass cont panel, ASN-75	X		
HSI	X		
FLIGHT INSTRUMENTS	No EMI Affect	EMI Affected Gnd Flt	Explanation
Radar altimeter	X		
Stabilator pos indicator	X		
VSI	X		
CIS mode select	X		
SAS 1	X		
SAS 2	X		
FPS	X		
Trim	X		
Go-around enable	X		
Cyclic trim release	X		
Cyclic stick trim	X		
ALR encoder	X		

FLIGHT INSTRUMENTS (CONT)	No EMI Affect	EMI Affected Gnd Flt	Explanation
HSI/VSI mode select (PLT)			
DPLR	X		
VOR/ILS	X		
BACK CRS	X		
FM HOME	X		
TURN RATE	X		
CRS HDG	X		
VERT GYRO	X		
BRG 2	X		
HSI/VSI Mode Select (COPLT)			
DPLR	X		
VOR/ILS	X		
BACK CRS	X		
FM HOME	X		
TURN RATE	X		
CRS HDG	X		
VERT GYRO	X		
BRG 2	X		
MISCELLANEOUS EQUIPMENT	No EMI Affect	EMI Affected Gnd Flt	Explanation
Blade deice	Not tested		Ambient tempera- ture was out of test lim- its.
Windshield anti-ice	X		
Pitot heat	X		
Vent blower	X		
Windshield wiper	X		
Heater	X		
APU	X		
Generator #1	X		
Generator #2	X		
Generator APU	X		
Air source heat start	X		
Tail wheel lock	X		
Gyro erect	X		

LIGHTING	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
Cockpit utility	X			
Cockpit flood	X			
Cabin dome	X			
Search light	X			
Search light control	X			
Landing light	X			
Flt instr lights (PLT)	X			
Flt instr lights (COPLT)	X			
Nonflight instr lights	X			
Console lights, upper	X			
Console lights, lower	X			
Position lights	X			
Formation lights	X			
Anticollision lights	X			
NVG lighting	X			

3.2.5 Battery life evaluation

Battery Life Evaluation Report Form

Nomenclature: Transport Incubator
Manufacturer: Ohio Medical Products
Model number: Air-Vac
Serial number: AKHD00465
Military item number: None

Options installed: None

Manufacturer battery life specification: Up to 3 hours while maintaining a temperature of 90°F in an ambient temperature of 70°F. Up to 1 hour in an ambient temperature of 40°F.

Overall performance: Not evaluated

Performance: The battery on the unit secured for laboratory testing would not charge. The unit used for flight testing provided 45 minutes operation while set at 90°F (58°F ambient).

Comments: None

3.2.6 Electrical safety test

Electrical Safety Test Report Form

Nomenclature: Transport Incubator
Manufacturer: Ohio Medical Products
Model number: Ohio Medical Air-Vac
Serial number: AKHD00465
Military item number: None

Options installed: None

Date of test: 20 Aug 91

Performance:

Grounding conductor resistance (milliohms): 50.7

Leakage current - Case to ground (microamperes):

unit off, grounded, normal polarity	0.0
unit off, ungrounded, normal polarity	7.2
unit off, ungrounded, reverse polarity	8.9
unit on, grounded, normal polarity	0.1
unit on, ungrounded, normal polarity	12.8
unit on, ungrounded, reverse polarity	12.6

MAXIMUM LIMITS:

ground resistance (milliohms):	150
current (microamperes)	
current (grounded, type A unit):	10
current (ungrounded, type A unit):	100
current (grounded, type B unit):	50
current (ungrounded, type B unit):	500

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.7 Human factors evaluation

Human Factors Evaluation Report Form

Nomenclature: Transport Incubator
Manufacturer: Ohio Medical Products
Model number: Air-Vac
Serial number: AKHD00465
Military item number: None

Options installed: None

Date of test: 11 Sep 91

Item configuration during test: Item prepared for operation.

Checklist for HFE

RESULTS

VISUAL DISPLAYS:

Satisfactory

- display type, format, content
- location of displays
- indicator lights
- scalar displays
- color coding
- legends and labels
- cathode ray tubes
- counters
- flags, go-no-go, center-null indicators

Comments: Power on is indicated by a red light. Red lighting is usually reserved for an error condition.

CONTROLS:

Satisfactory

- location
- characteristics of controls
- labeling
- control - display relationships

Comments: Temperature control not graduated in degrees.

TIME REQUIRED TO PREPARE FOR OPERATION (list in comment)

Comments: approximately 2 minutes.

MAINTAINABILITY:

Satisfactory

- component location
- component characteristics
- rests and stands
- covers, cases, access doors
- handles
- lubrication
- component mounting
- cord storage provisions
- external accessibility
- internal accessibility
- list special tools required
- list realistic inspection requirements
- list realistic inspection intervals

Comments: Operational checks should be performed as used and maintenance inspection every 6 months.

CONDUCTORS:

Satisfactory

- binding and securing
- length
- protection
- routing
- conductor coding
- fabrication
- connectors

Comments: Unit received for evaluation had cracked insulation on the power cord. Another (newer) unit was checked and this was not present.

FASTENERS:

Satisfactory

- access through inspection panel covers
- enclosure fasteners
- device mounting bolts and fasteners

Comments: Several access door latches were broken on an older unit, but a new unit had all fasteners intact.

TEST POINTS:

Satisfactory

general
location and mounting
test point labeling and coding

Comments: None

TEST EQUIPMENT:

Satisfactory

general
equipment self-test
indicators (list in comments)
controls
positive indication of proper operation

Comments: None

FUSES AND CIRCUIT BREAKERS:

Satisfactory

external accessibility
easy replacement or reset by operator

Comments: None

LABELS AND CODING:

Satisfactory

placed above controls and displays
near or on the items they identify
not obscured by other equipment components
describe the function of the items they identify
readable from normal operating distance
conspicuous placards adjacent to hazardous items

Comments: None

SAFETY:

Satisfactory

manual
materials
fire and explosive protection
operator protection from mechanical hazards
patient protection from mechanical hazards
electrical safety (operator and patient)

Comments: The first unit tested had a different
battery model than the one described
in the manual and the power cord was
damaged. This was repaired prior to use.

3.2.8 Altitude test

Altitude Test Report Form

Nomenclature: Transport Incubator
Manufacturer: Ohio Medical Products
Model number: Ohio Medical Air-Vac
Serial number: AKHD00465
Military item number: None

Options installed: None

Date of test: 28 Aug 91

Item configuration during test: Item sitting on chamber floor.

Performance test criteria: Accurate maintenance of selected temperature.

Ambient conditions outside chamber:

Temperature	20°C
Humidity	87% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	None

IN-TEST DATA

Time of test start: 0800

POSTTEST DATA

Posttest performance check (complete check of item and accessories):

Time of test end : 0910

Item functional (based on performance test criteria): Yes

Deviation from pretest : None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.9 Vibration test

Vibration Test Report Form

Nomenclature: Transport Incubator
Manufacturer: Ohio Medical Products
Model number: Ohio Medical Air-Vac
Serial number: AKHD00465
Military item number: None

Options installed: None

Date of test: 4 Sep 91

Item configuration during test: Item strapped down on vibration table fixture; ac operation.

Performance test criteria: Accurate maintenance of selected temperature.

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	None

Ambient conditions

Temperature	19°C
Humidity	75% RH
Barometric pressure	1 atm

IN-TEST DATA

Data and performance checks during test:

Time at first check:

X: 1235 Y: 1346 Z: 0850

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Time at second check:

X: 1335

Y: 1443

Z: 0945

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

POSTTEST DATA

Time at test end:

X: 1336

Y: 1447

Z: 0950

Posttest performance check (complete check of item and accessories):

Item functional (based on performance test criteria): Yes

Item intact: Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: Test times for the three axes are on different days.

3.2.10 High temperature test

High Temperature Test (Equipment Operating) Report Form

Nomenclature: Transport Incubator
Manufacturer: Ohio Medical Products
Model number: Ohio Medical Air-Vac
Serial number: AKHD00465
Military item number: None

Options installed: None

Date of test: 29 Aug 91

Item configuration during test: Unit was sitting on chamber floor, operating on ac power.

Performance test criteria: Accurate delivery of correct volume of fluid at a rate of 60 mL/h (measured).

Ambient conditions outside chamber:

Temperature	25°C
Humidity	54% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check :

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	None
distance from north wall (meters)	0.504
distance from south wall (meters)	1.080
distance from east wall (meters)	1.07
distance from west wall (meters)	1.17
distance from ceiling (meters)	0.97
distance from floor (meters)	0.79

IN-TEST DATA

Time of test start: 0835

Performance checks during test:

First check:

Time: 0905
Temperature: 49°C
Humidity: 15% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria):
Yes, all ok
Deviation from pretest: Incubator control set for 90°F,
canopy temperature 107°F and overtemp alarm not
activated.

Second check:

Time: 0935
Temperature: 49°C
Humidity: 15% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria):
Yes, all ok
Deviation from pretest: Canopy temperature 114°F and
no overtemp alarm indication.

Third check:

Time: 1005
Temperature: 49°C
Humidity: 15% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria):
Yes, all ok
Deviation from pretest: Canopy temperature 116°F and
no overtemp alarm indication.

POSTTEST DATA

Posttest performance check:
(complete check of item and accessories)

Time of test end: 1100
Item functional (based on performance test criteria):
Yes, all ok
Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: High temperature alarm activated at
120°F indicating improper calibration on this unit.

3.2.11 High temperature storage test

High Temperature Test (Equipment in Storage) Report Form

Nomenclature: Transport Incubator
Manufacturer: Ohio Medical Products
Model number: Ohio Medical Air-Vac
Serial number: AKHD00465
Military item number: None

Options installed: None

Date of test: 5 Sep 91

Item configuration during test: Sitting on chamber floor, in storage, not operating.

Performance test criteria: Consistent and accurate maintenance of selected incubator temperature.

Ambient conditions outside chamber:

Temperature	24C
Humidity	55% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	None
distance from north wall (meters)	0.56
distance from south wall (meters)	0.56
distance from east wall (meters)	1.07
distance from west wall (meters)	1.17
distance from ceiling (meters)	0.97
distance from floor (meters)	0.79

Time of test start: 0930

POSTTEST DATA

Posttest performance check:
(complete check of item and accessories)

Time of test end: 1230
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Comments on item setup or checks:
The unit was allowed to cool for 1 hour at ambient conditions before the posttest performance check was completed.

Comments on test run (including interruptions): None

Comments on other data: None

3.2.12 Low temperature test

Low Temperature Test (Equipment Operating) Report Form

Nomenclature: Transport Incubator
Manufacturer: Ohio Medical Products
Model number: Ohio Medical Air-Vac
Serial number: AKHD00465
Military item number: None

Options installed: None

Date of test: 29 Aug 91

Item configuration during test: Sitting on chamber floor.

Performance test criteria: Accurate maintenance of selected incubator temperature.

Ambient conditions outside chamber:

Temperature	23°C
Humidity	45% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Pass

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	None
distance from north wall (meters)	0.56
distance from south wall (meters)	0.56
distance from east wall (meters)	1.07
distance from west wall (meters)	1.17
distance from ceiling (meters)	0.97
distance from floor (meters)	0.79

Time of test start: 1200

Performance checks during test:

First check:

Time: 1230
Temperature: 0°C
Humidity: NA
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: Selected temperature is 90°F,
incubator temperature is 77°F.

Second check:

Time: 1300
Temperature: 0°C
Humidity: NA
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: Selected temperature is 90°F,
incubator temperature is 75°F.

Third check:

Time: 1330
Temperature: 0°C
Humidity: NA
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: Selected temperature is 90°F,
incubator temperature is 73°F.

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1409
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Comments on item setup or checks: Incubator temperature was
90°F at beginning of test.

Comments on test run (including interruptions): Incubator
unable to maintain selected temperature.

Comments on other data: None

3.2.13 Low temperature storage test

Low Temperature Test (Equipment in Storage) Report Form

Nomenclature: Transport Incubator
Manufacturer: Ohio Medical Products
Model number: Ohio Medical Air-Vac
Serial number: AKHD00465
Military item number: None

Options installed: None

Date of test: 6 Sep 91

Item configuration during test: Sitting on chamber floor, not operating, in storage.

Performance test criteria: Accurate maintenance of selected incubator temperature.

Ambient conditions outside chamber:

Temperature	25°C
Humidity	55% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	None
distance from north wall (meters)	0.56
distance from south wall (meters)	0.56
distance from east wall (meters)	1.07
distance from west wall (meters)	1.17
distance from ceiling (meters)	0.97
distance from floor (meters)	0.79

Time of test start:	0900
Midtest time:	1200
Midtest temperature:	-46°C

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1545

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.14 Humidity test

Humidity Test Report Form

Nomenclature: Transport Incubator
Manufacturer: Ohio Medical Products
Model number: Model Air-Vac
Serial number: AKHD00465
Military item number: None

Options installed: None

Date of test: 30 Aug 91

Item configuration during test: The unit was sitting on the chamber floor, operating on ac power.

Performance test criteria: Accurate maintenance of selected temperature.

Ambient conditions outside chamber:

Temperature	25°C
Humidity	57% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	None
distance from north wall (meters)	0.56
distance from south wall (meters)	0.56
distance from east wall (meters)	1.07
distance from west wall (meters)	1.17
distance from ceiling (meters)	0.97
distance from floor (meters)	0.79

IN-TEST DATA

Time of test start: 0827

Performance checks during test:

First check:

Time: 0915
Temperature: 29.5°C
Humidity: 95% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Second check:

Time: 1000
Temperature: 29.5°C
Humidity: 95% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Third check:

Time: 1045
Temperature: 29.5°C
Humidity: 95% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Fourth check:

Time: 1130
Temperature: 29.5°C
Humidity: 95% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Fifth check:

Time: 1215
Temperature: 29.5°C
Humidity: 95% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1315

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.15 Electromagnetic characteristics test

Electromagnetic Characteristics Testing Evaluation of Performance

T & E Item Number: 33

Date: 27 Aug 91

Nomenclature: Transport Incubator
Manufacturer: Ohio Medical Products
Model number: Air-Vac
Serial number: AKHD00465
Military item number: NA

Conducted Emissions Tests

CE01 Testing configuration(s): NA
 Performance (pass/fail): NA

Comments: No dc conductors

CE02 Testing configuration(s): Operating on copper
 work bench.
 Performance (pass/fail): Pass

Comments: No signal failures.

CE04 Testing configuration(s): Operating on copper
 work bench.
 Performance (pass/fail): Pass

Comments: No signal failures.

Conducted Susceptibility Tests

CS02 Testing configuration(s): Operating on test
 bench, connected to test jig.
 Performance (pass/fail): Pass

Comments: Not susceptible to test signals on
power conductors.

CS06 Testing configuration(s): Operating on counter
 top.
 Performance (pass/fail): Pass

 Comments: Not susceptible to test spikes

Radiated Emissions Tests

RE02 Testing configuration(s): Operating on wooden
 test stand in the EMC chamber, ac power.
 Performance (pass/fail): Fail

 Comments: BB failures 0.4 to 1.7 dB over specifi-
 cations in range 150 to 245 KHz; 0.8 to 9.5 dB
 over specifications in range 4.7 to 17.97 MHz;
 and 1.7 to 13.4 dB over specifications in range
 23.44 to 27.81 MHz.

Radiated Susceptibility Tests

RS03 Testing configuration(s): Operating on the wooden
 test stand in the EMC chamber.
 Performance (pass/fail): Pass

 Comments: Battery power tests terminated at 200
 MHz because battery could not be recharged.
 Remainder of test on ac power. Not susceptible
 to test signals.

3.3 CRITERIA, SIGNIFICANT PROBLEMS, AND SUGGESTED IMPROVEMENTS

3.3.1 Criteria

Item			<u>Applicable</u>
<u>No.</u>	<u>Criteria (source)</u>	<u>Remarks</u>	<u>subparagraph</u>
1	The physical inventory is conducted solely for investigation and documentation.	NA	2.1.2.1
2	The Ohio Medical Air-Vac will display consistent and accurate performance.	partially met	2.1.2.2
3	Verify manufacturer's specified full power internal battery life expectancy of 10 hours.	not evaluated	2.2.2
4	The Ohio Medical Air-Vac will meet the limits established in NFPA 99 for electrical safety of medical equipment.	met	2.3.2
5	The Ohio Medical Air-Vac will be rated satisfactory in all major categories of the evaluation. These include: Visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.	met	2.4.2
6	The Ohio Medical Air-Vac will demonstrate proper operation - while exposed to an altitude equivalency of 15,000 feet above sea level.	met	2.5.2
7	The Ohio Medical Air-Vac will remain operational while exposed to vibrational stresses.	met	2.6.2
8	The Ohio Medical Air-Vac will remain operational during the high temperature operation check.	partially met	2.7.2.1

9	The Ohio Medical Air-Vac will remain operational after the high temperature storage.	met	2.7.2.2
10	The Ohio Medical Air-Vac will remain operational during the low temperature operation check.	partially met	2.8.2.1
11	The Ohio Medical Air-Vac will remain operational after the low temperature storage.	met	2.8.2.2
12	The Ohio Medical Air-Vac will remain operational while exposed to a high humidity.	met	2.9.2
13	The Ohio Medical Air-Vac will not produce emissions in excess of the limits set forth in MIL-STD-461A Notice 4, paragraph 6.13.	partially met	2.10.2.1
14	The Ohio Medical Air-Vac will not malfunction when it is subjected to radiated fields as specified in MIL-STD-461A, Notice 4, paragraph 6.20.	met	2.10.2.2
15	The Ohio Medical Air-Vac will not conduct emissions in excess of the limits set forth in MIL-STD-461A, Notice 4, paragraph 6.2.	met	2.10.2.3
16	The Ohio Medical Air-Vac will not malfunction when it is subjected to conducted emissions as specified in MIL-STD-461A, Notice 4, paragraphs 6.7 and 6.10.	met	2.10.2.4
17	The flight surgeon will be able to operate the Ohio Medical Air-Vac without physical or functional restrictions aboard the aircraft.	met	2.11.2.1

- | | | | |
|----|--|-----|----------|
| 18 | The Ohio Medical Air-Vac will not radiate EMI to disrupt or interfere with the other equipment or systems aboard the aircraft. | met | 2.12.2.2 |
| 19 | The aircraft will not radiate EMI to disrupt or interfere with the Ohio Medical Air-Vac. | met | 2.12.2.3 |

3.3.2 Significant problems which require corrective action

None

3.3.3 Suggested improvements

None

3.4 REFERENCES

3.4.1 Department of Defense. 1971. EMI characteristics, requirements for equipment. Washington, DC. MIL-STD-461A, Notice 4. February.

3.4.2 Department of Defense. 1971. EMI characteristics, measurement of. Washington, DC. MIL-STD-462, Notice 3. February.

3.4.3 Department of Defense. 1983. Environmental test methods and engineering guidelines. Washington, DC. MIL-STD-810D. July.

3.4.4 Department of the Army. 1987. Maintenance management procedures for medical equipment. Washington, DC. TB 38-750-2. April.

3.4.5 Underwriters Laboratory's, Inc. 1978. Standard for safety, medical and dental equipment. Chicago, Illinois. UL-544.

3.4.6 Department of Defense. 1989. Human engineering design criteria for military systems, equipment, and facilities. Washington, DC. MIL-STD-1472D. March.

3.4.7 Association for the Advancement of Medical Instruments. 1988. Human factors engineering guidelines and preferred practices for the design of medical devices. Arlington, Virginia. AAMI-HE-1988. February.

3.4.8 National Fire Protection Association. 1987. Standard for health care facilities. Quincy, Massachusetts. NFPA 99. February.

3.4.9 Department of the Army. 1982. Environmental protection and enhancement. Washington, DC. AR 200-1. June.

3.5 ABBREVIATIONS

ac	alternate current
AVSCOM	Army Aviation Systems Command
AWR	airworthiness release
BB	broadband
CAAF	Cairns Army Airfield
dc	direct current
EMC	electromagnetic compatibility
EMI	electromagnetic interference
fpm	feet per minute
GFE	government furnished equipment
Gpk	gravity, peak
G(rms)	gravity (root mean square)
Hz	hertz
IAW	in accordance with
ITOP	in-flight test operating procedure
IV	intravenous
kHz	kilohertz
LCD	liquid crystal display
LED	light emitting diode
LISN	line impedance stabilization network
MEDEVAC	medical evacuation
MHz	megahertz
MIL-STD	military standard
mL	milliliter
mm	millimeter
mmHg	millimeters of Mercury
MSL	mean sea level
NFPA	National Fire Prevention Association
NB	narrowband
NBC	nuclear, biological and chemical
NOE	nap-of-the-earth
NVG	night vision goggle
RF	radio frequency
RFI	radio frequency interference
RH	relative humidity

TB	technical bulletin
TFT	technical feasibility testing
T & E	test and evaluation
UES	Universal Energy Systems, Inc.
USAARL	U.S. Army Aeromedical Research Laboratory
V/m	volts per meter

3.6 LIST OF MANUFACTURERS

- 3.6.1 Ohio Medical Products
4565 Winters Chapel Road
Atlanta, GA 30340
- 3.6.2 Neurodyne-Dempsey, Inc.
200 Arrowhead Drive
Carson City, NV 89701
- 3.6.3 Tenney Engineering, Inc.
1090 Springfield Road
P.O. box 3142
Union, NJ 07083
- 3.6.4 Unholtz-Dickey Corporation
6 Brookside Drive
Wallingford, CT 06492
- 3.6.5 Solar Electronics Company
901 North Highland Avenue
Hollywood, CA 90038
- 3.6.6 Tektronix, Inc.
P.O. Box 500
Beaverton, OR 97077

3.7 DISTRIBUTION LIST

Commander, U.S. Army Natick Research,
Development and Evaluation Center
ATTN: STRNC-MIL (Documents
Librarian)
Natick, MA 01760-5040

Commander
U.S. Army Aviation Systems Command
ATTN: AMSAV-ECU
4300 Goodfellow Bouvelard
St. Louis, MO 63120-1790

Commander/Director
U.S. Army Combat Surveillance
and Target Acquisition Lab
ATTN: DELCS-D
Fort Monmouth, NJ 07703-5304

Commander
10th Medical Laboratory
ATTN: Audiologist
APO New York 09180

Naval Air Development Center
Technical Information Division
Technical Support Detachment
Warminster, PA 18974

Commanding Officer, Naval Medical
Research and Development Command
National Naval Medical Center
Bethesda, MD 20814-5044

Deputy Director, Defense Research
and Engineering
ATTN: Military Assistant
for Medical and Life Sciences
Washington, DC 20301-3080

Commander, U.S. Army Research
Institute of Environmental Medicine
Natick, MA 01760

U.S. Army Avionics Research
and Development Activity
ATTN: SAVAA-P-TP
Fort Monmouth, NJ 07703-5401

U.S. Army Communications-Electronics
Command
ATTN: AMSEL-RD-ESA-D
Fort Monmouth, NJ 07703

Library
Naval Submarine Medical Research Lab
Box 900, Naval Sub Base
Groton, CT 06349-5900

Commander
Man-Machine Integration System
Code 602
Naval Air Development Center
Warminster, PA 18974

Commander
Naval Air Development Center
ATTN: Code 602-B (Mr. Brindle)
Warminster, PA 18974

Commanding Officer
Harry G. Armstrong Aerospace
Medical Research Laboratory
Wright-Patterson
Air Force Base, OH 45433

Director
Army Audiology and Speech Center
Walter Reed Army Medical Center
Washington, DC 20307-5001

Commander, U.S. Army Institute
of Dental Research
ATTN: Jean A. Setterstrom, Ph. D.
Walter Reed Army Medical Center
Washington, DC 20307-5300

Naval Air Systems Command
Technical Air Library 950D
Room 278, Jefferson Plaza II
Department of the Navy
Washington, DC 20361

Naval Research Laboratory Library
Shock and Vibration
Information Center, Code 5804
Washington, DC 20375

Director, U.S. Army Human
Engineering Laboratory
ATTN: Technical Library
Aberdeen Proving Ground, MD 21005

Commander, U.S. Army Test
and Evaluation Command
ATTN: AMSTE-AD-H
Aberdeen Proving Ground, MD 21005

Director
U.S. Army Ballistic
Research Laboratory
ATTN: DRXBR-OD-ST Tech Reports
Aberdeen Proving Ground, MD 21005

Commander
U.S. Army Medical Research
Institute of Chemical Defense
ATTN: SGRD-UV-AO
Aberdeen Proving Ground,
MD 21010-5425

Commander U.S. Army Medical
Research and Development Command
ATTN: SGRD-RMS (Ms. Madigan)
Fort Detrick, Frederick, MD 21702-5012

Director
Walter Reed Army Institute of Research
Washington, DC 20307-5100

HQ DA (DASG-PSP-O)
5109 Leesburg Pike
Falls Church, VA 22041-3258

Naval Research Laboratory
Library Code 1433
Washington, DC 20375

Harry Diamond Laboratories
ATTN: Technical Information Branch
2800 Powder Mill Road
Adelphi, MD 20783-1197

U.S. Army Materiel Systems
Analysis Agency
ATTN: AMXSYP-PA (Reports Processing)
Aberdeen Proving Ground
MD 21005-5071

U.S. Army Ordnance Center
and School Library
Simpson Hall, Building 3071
Aberdeen Proving Ground, MD 21005

U.S. Army Environmental
Hygiene Agency
Building E2100
Aberdeen Proving Ground, MD 21010

Technical Library Chemical Research
and Development Center
Aberdeen Proving Ground, MD
21010-5423

Commander
U.S. Army Medical Research
Institute of Infectious Disease
SGRD-UIZ-C
Fort Detrick, Frederick, MD 21702

Director, Biological
Sciences Division
Office of Naval Research
600 North Quincy Street
Arlington, VA 22217

Commander
U.S. Army Materiel Command
ATTN: AMCDE-XS
5001 Eisenhower Avenue
Alexandria, VA 22333

Commandant
U.S. Army Aviation
Logistics School ATTN: ATSQ-TDN
Fort Eustis, VA 23604

Headquarters (ATMD)
U.S. Army Training
and Doctrine Command
Fort Monroe, VA 23651

Structures Laboratory Library
USARTL-AVSCOM
NASA Langley Research Center
Mail Stop 266
Hampton, VA 23665

Naval Aerospace Medical
Institute Library
Building 1953, Code 03L
Pensacola, FL 32508-5600

Command Surgeon
HQ USCENTCOM (CCSG)
U.S. Central Command
MacDill Air Force Base FL 33608

Air University Library
(AUL/LSE)
Maxwell Air Force Base, AL 36112

U.S. Air Force Institute
of Technology (AFIT/LDEE)
Building 640, Area B
Wright-Patterson
Air Force Base, OH 45433

Henry L. Taylor
Director, Institute of Aviation
University of Illinois-Willard Airport
Savoy, IL 61874

COL Craig L. Urbauer, Chief
Office of Army Surgeon General
National Guard Bureau
Washington, DC 50310-2500

Commander
U.S. Army Aviation Systems Command
ATTN: SGRD-UAX-AL (MAJ Gillette)
4300 Goodfellow Blvd., Building 105
St. Louis, MO 63120

U.S. Army Aviation Systems Command
Library and Information Center Branch
ATTN: AMSAV-DIL
4300 Goodfellow Boulevard
St. Louis, MO 63120

Federal Aviation Administration
Civil Aeromedical Institute
Library AAM-400A
P.O. Box 25082
Oklahoma City, OK 73125

Commander
U.S. Army Academy
of Health Sciences
ATTN: Library
Fort Sam Houston, TX 78234

Commander
U.S. Army Institute of Surgical Research
ATTN: SGRD-USM (Jan Duke)
Fort Sam Houston, TX 78234-6200

AAMRL/HEX
Wright-Patterson
Air Force Base, OH 45433

John A. Dellinger,
Southwest Research Institute
P. O. Box 28510
San Antonio, TX 78284

Product Manager
Aviation Life Support Equipment
ATTN: AMCPM-ALSE
4300 Goodfellow Boulevard
St. Louis, MO 63120-1798

Commander
U.S. Army Aviation
Systems Command
ATTN: AMSAV-ED
4300 Goodfellow Boulevard
St. Louis, MO 63120

Commanding Officer
Naval Biodynamics Laboratory
P.O. Box 24907
New Orleans, LA 70189-0407

Assistant Commandant
U.S. Army Field Artillery School
ATTN: Morris Swott Technical Library
Fort Sill, OK 73503-0312

Commander
U.S. Army Health Services Command
ATTN: HSOP-SO
Fort Sam Houston, TX 78234-6000

Director of Professional Services
HQ USAF/SGDT
Bolling Air Force Base, DC 20332-6188

U.S. Army Dugway Proving Ground
Technical Library, Building 5330
Dugway, UT 84022

U.S. Army Yuma Proving Ground
Technical Library
Yuma, AZ 85364

AFFTC Technical Library
6510 TW/TSTL
Edwards Air Force Base,
CA 93523-5000

Commander
Code 3431
Naval Weapons Center
China Lake, CA 93555

Aeromechanics Laboratory
U.S. Army Research and Technical Labs
Ames Research Center, M/S 215-1
Moffett Field, CA 94035

Sixth U.S. Army
ATTN: SMA
Presidio of San Francisco, CA 94129

Commander
U.S. Army Aeromedical Center
Fort Rucker, AL 36362

U.S. Air Force School
of Aerospace Medicine
Strughold Aeromedical Library Technical
Reports Section (TSKD)
Brooks Air Force Base, TX 78235-5301

Dr. Diane Damos
Department of Human Factors
ISSM, USC
Los Angeles, CA 90089-0021

U.S. Army White Sands
Missile Range
ATTN: STEWS-IM-ST
White Sands Missile Range, NM 88002

U.S. Army Aviation Engineering
Flight Activity
ATTN: SAVTE-M (Tech Lib) Stop 217
Edwards Air Force Base, CA 93523-5000

Ms. Sandra G. Hart
Ames Research Center
MS 262-3
Moffett Field, CA 94035

Commander, Letterman Army Institute
of Research
ATTN: Medical Research Library
Presidio of San Francisco, CA 94129

COL Eugene S. Channing, O.D.
Brooke Army Medical Center
ATTN: HSHE-EAH-O
Fort Sam Houston, TX 78234-6200

Commander
U.S. Army Medical Materiel
Development Activity
Fort Detrick, Frederick, MD 21702-5009

Commander
U.S. Army Aviation Center
Directorate of Combat Developments
Building 507
Fort Rucker, AL 36362

U. S. Army Research Institute
Aviation R&D Activity
ATTN: PERI-IR
Fort Rucker, AL 36362

Commander
U.S. Army Safety Center
Fort Rucker, AL 36362

U.S. Army Aircraft Development
Test Activity
ATTN: STEBG-MP-P
Cairns Army Air Field
Fort Rucker, AL 36362

Commander U.S. Army Medical Research
and Development Command
ATTN: SGRD-PLC (COL Sedge)
Fort Detrick, Frederick, MD 21702

MAJ John Wilson
TRADOC Aviation LO
Embassy of the United States
APO New York 09777

Netherlands Army Liaison Office
Building 602
Fort Rucker, AL 36362

British Army Liaison Office
Building 602
Fort Rucker, AL 36362

Italian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Directorate of Training Development
Building 502
Fort Rucker, AL 36362

Chief
USAHEL/USAAVNC Field Office
P. O. Box 716
Fort Rucker, AL 36362-5349

Commander U.S. Army Aviation Center
and Fort Rucker
ATTN: ATZQ-CG
Fort Rucker, AL 36362

Commander/President
TEXCOM Aviation Board
Cairns Army Air Field
Fort Rucker, AL 36362

MAJ Terry Newman
Canadian Army Liaison Office
Building 602
Fort Rucker, AL 36362

German Army Liaison Office
Building 602
Fort Rucker, AL 36362

LTC Patrice Cottebrune
French Army Liaison Office
USAAVNC (Building 602)
Fort Rucker, AL 36362-5021

Brazilian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Australian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Dr. Garrison Rapmund
6 Burning Tree Court
Bethesda, MD 20817

Commandant Royal Air Force
Institute of Aviation Medicine
Farnborough Hants UK GU14 6SZ
Dr. A. Kornfield, President
Biosearch Company
3016 Revere Road
Drexel Hill, PA 29026

Commander
U.S. Army Biomedical Research
and Development Laboratory
ATTN: SGRD-UBZ-I
Fort Detrick, Frederick, MD 21702

Defense Technical Information Center
Cameron Station
Alexandra, VA 22313

Commander, U.S. Army Foreign Science
and Technology Center
AIFRTA (Davis)
220 7th Street, NE
Charlottesville, VA 22901-5396

Director,
Applied Technology Laboratory
USARTL-AVSCOM
ATTN: Library, Building 401
Fort Eustis, VA 23604

U.S. Army Training
and Doctrine Command
ATTN: Surgeon
Fort Monroe, VA 23651-5000

Aviation Medicine Clinic
TMC #22, SAAF
Fort Bragg, NC 28305

U.S. Air Force Armament
Development and Test Center
Eglin Air Force Base, FL 32542

Commander, U.S. Army Missile
Command
Redstone Scientific Information Center
ATTN: AMSMI-RD-CS-R/ILL
Documents Redstone Arsenal, AL 35898

U.S. Army Research and Technology
Laboratories (AVSCOM)
Propulsion Laboratory MS 302-2
NASA Lewis Research Center
Cleveland, OH 44135

Dr. H. Dix Christensen
Bio-Medical Science Building, Room 753
Post Office Box 26901
Oklahoma City, OK 73190

Dr. Christine Schlichting
Behavioral Sciences Department
Box 900, NAVUBASE NLON
Groton, CT 06349-5900

Commandant
Academy of Health Sciences
ATTN: HSHA-COM (LTC Huether)
Fort Sam Houston, TX 78234

U.S. Air Force Armament
Development and Test Center
Eglin Air Force Base, FL 32542

COL Eugene S. Channing, O.D.
Brooke Army Medical Center
ATTN: HSHE-EAH-O
Fort Sam Houston, TX 78234-6200